


GC CAPILLARY COLUMN APPLICATIONS: TRIGLYCERIDES

Butter Triglycerides

COLUMN: 007-65HT (65% Phenyl)
Methylpolysiloxane
15M. x 0.25mm I.D. x
0.1µm film
Cat. No.: 400-65HT-15-0.1F
Temperature: 250° (5°/min.) - 365°C
Injector: 370°C
Detector: 400°C, FID
Carrier Gas: 35 cm/sec., Helium

1. cholesterol
2. ppp
3. pps
4. ppo
5. pss
6. pso
7. poo
8. ssp
9. soo
10. ooo

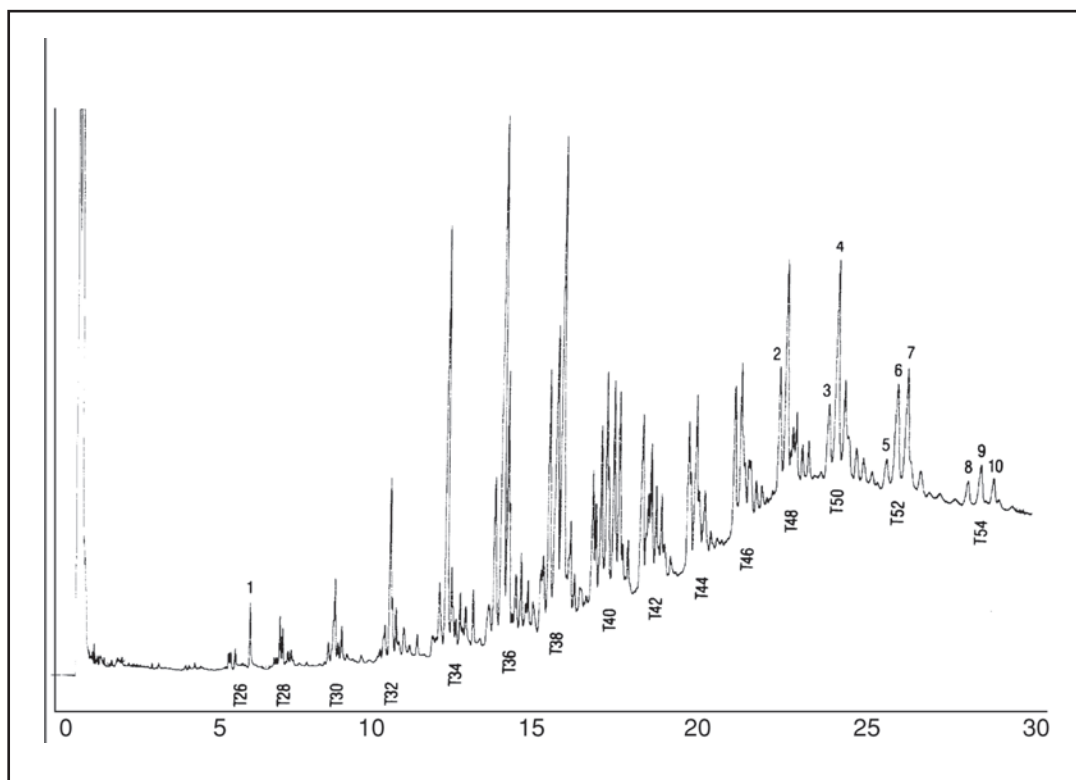

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com