


GC CAPILLARY COLUMN APPLICATIONS: FLAVORS/FRAGRANCES

Celery Seed Oil

COLUMN: BTR-CW, Carbowax/
BTR -Oxygen Resistant
30M. x 0.25mm I.D. x
0.25µm film
Cat. No.: BTR-CW-30-0.25F
Temperature: 60° (5°/min.) - 240°C
Injector: 225°C
Detector: 250°C, FID
Carrier Gas: 27 cm/sec., Helium

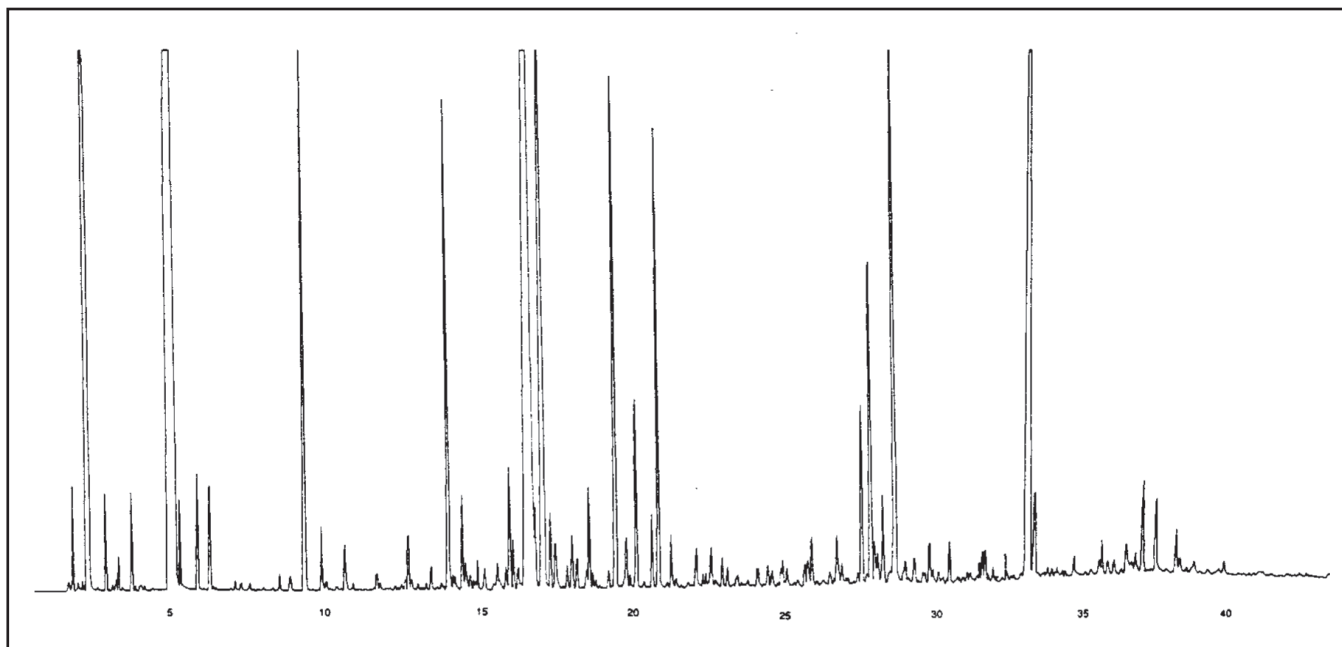

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com