


GC CAPILLARY COLUMN APPLICATIONS: TRIGLYCERIDES

Cocoa Butter

COLUMN: 007-65HT (65% Phenyl)
Methylpolysiloxane
15M. x 0.25mm I.D. x
0.1µm film
Cat. No.: 400-65HT-15-0.1F
Temperature: 350° (0.5min. hold)
(1°/min.) - 365°C
Injector: 400°C
Detector: 400°C, FID
Carrier Gas: 35cm/sec., Helium

1. ppp
2. mop
3. pps
4. pop
5. plp
6. pss
7. pos
8. poo
9. pls
10. sss
11. sos
12. soo
13. soc

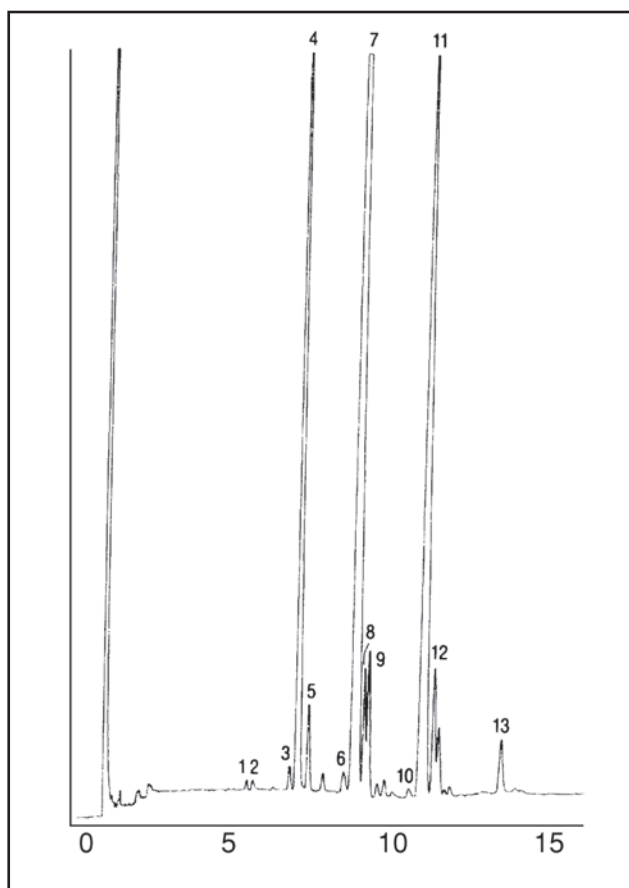

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com