

GC CAPILLARY COLUMN APPLICATIONS: TRIGLYCERIDES

Corn Oil

COLUMN: 400-65HT, (65% Phenyl)
Methylpolysiloxane -
Aluminum-Clad
15M. x 0.25mm I.D. x 0.1µm film

Cat. No.: 400-65HT-15-0.1F


Temperature: 350° (0.5 min. hold)
(1°/min.) - 365°C

Injector: 400°C

Detector: 400°C, FID

Carrier Gas: 35 cm/sec., Helium

1. pop
2. plp
3. pos
4. poo
5. pls
6. plo
7. pll
8. sos
9. soo
10. ooo
11. slo
12. olo
13. sll
14. oll
15. ll
16. ll

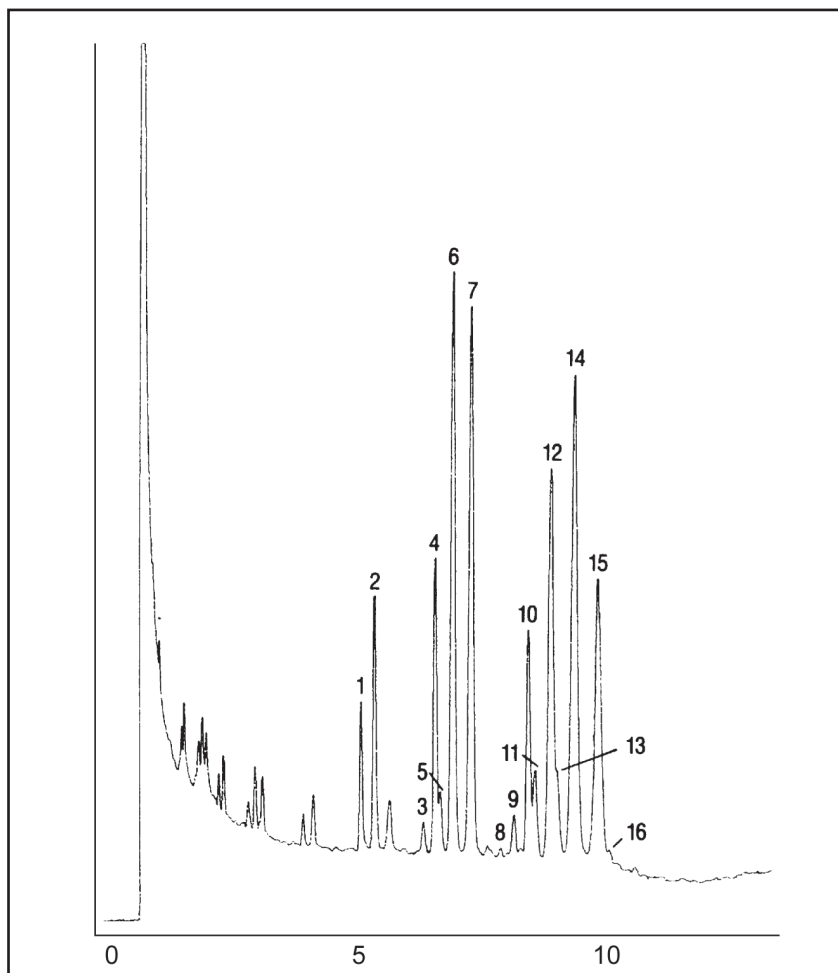

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com