


GC CAPILLARY COLUMN APPLICATIONS: ENVIRONMENTAL

DDT Isomers

COLUMN: 007-608,
Cyanopropyl Phenyl
Methylpolysiloxane
30M. x 0.53mm I.D. x
0.8 μ m film
Cat. No.: 007-608-30V-0.8F
Temperature: 140° (6°/min.) -260°C
Injector: 240°C
Detector: 300°C , ECD
Carrier Gas: 5 ml/min., Helium

1. o,p'-DDE
2. p,p'-DDE
3. o,p'-DDD
4. m,p'-DDD
5. p,p'-DDD
6. o,p'-DDT
7. p,p'-DDT

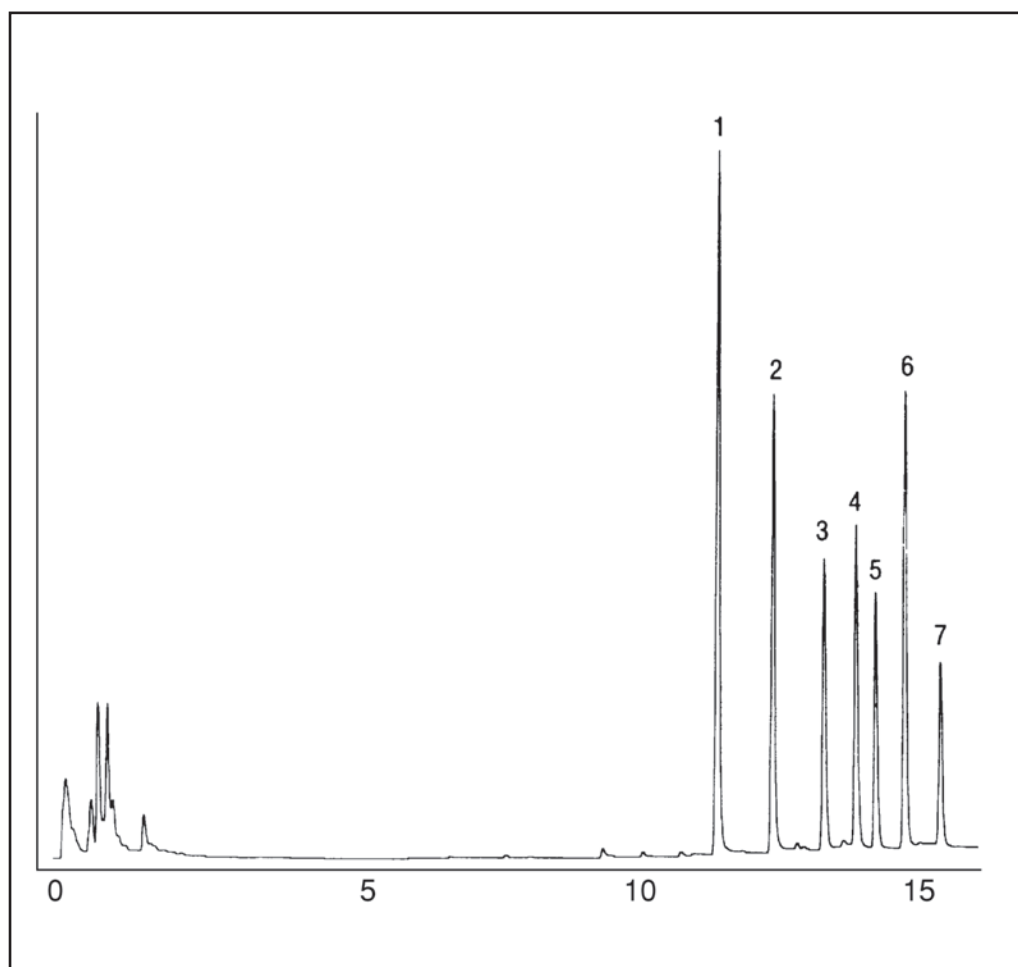

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com