


GC CAPILLARY COLUMN APPLICATIONS: ULTRA ALLOY™ STAINLESS STEEL

Solvent Mixture

COLUMN: UAC-1,
Dimethylpolysiloxane
60M. x 0.25mm I.D. x
1.0µm film
Cat. No.: UAC-1-60-1.0F
Temperature: 50° (15 min. hold)
(10°/min.) - 250°C
Injector: 200°C
Detector: 300°C

1. methanol
2. ethanol
3. acetone
4. isopropyl alcohol
5. dichloromethane
6. ethylacetate
7. chloroform
8. tetrahydrofuran
9. n-butyl alcohol
10. pyridine
11. dimethylformamide
12. cyclohexanone
13. phenethylamine

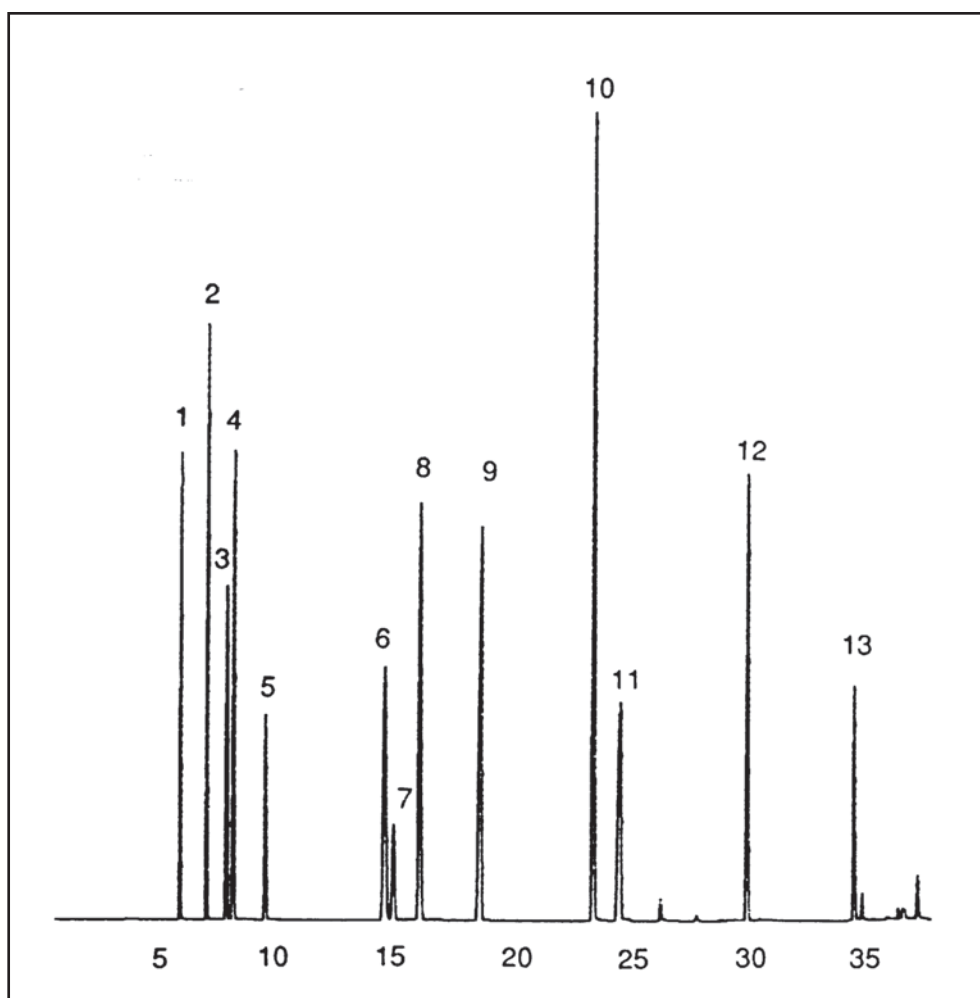

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com