


GC CAPILLARY COLUMN APPLICATIONS: ULTRA ALLOY™ STAINLESS STEEL

Amines

COLUMN: UAC-1,
Dimethylpolysiloxane
10M. x 0.25mm I.D. x 0.25µm film
Cat. No.: UAC-1-10-0.25F
Temperature: 100° (1 min. hold)
(20°/min.) - 300°C
Injector: 300°C
Detector: 400°C, FID
Carrier Gas: 30 cm/sec, Helium

1. di - n-butylamine
2. 4-dimethylaminopyridine
3. m-toluenediamine

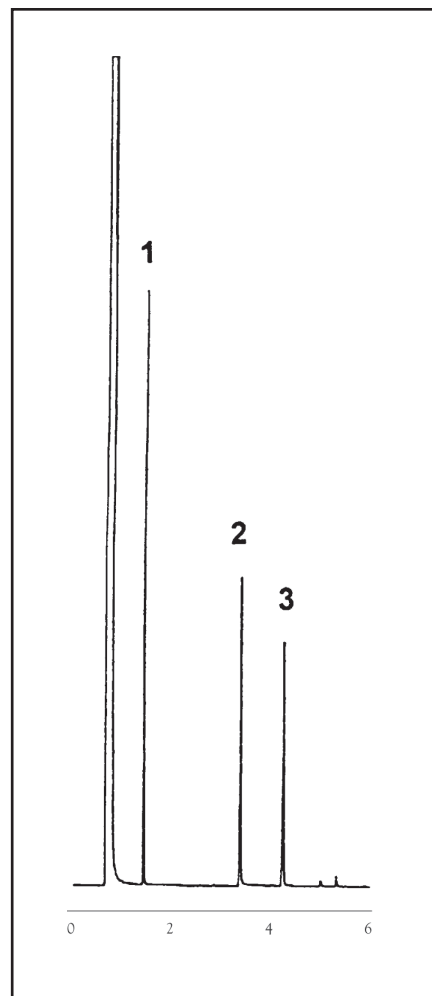

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com