


# GC CAPILLARY COLUMN APPLICATIONS: ULTRA ALLOY™ STAINLESS STEEL

## BTEX

COLUMN: UAC-CW, Carbowax  
50M. x 0.50mm I.D. x 1.0µm film  
Cat. No.: UAC-CW-50V-1.0F  
Temperature: 60°C  
Carrier Gas: H2  
Split Vent Flow: 100 ml/min.  
Injection: 0.1µL

1. benzene
2. toluene
3. ethylbenzene
4. p-xylene
5. m-xylene
6. o-xylene

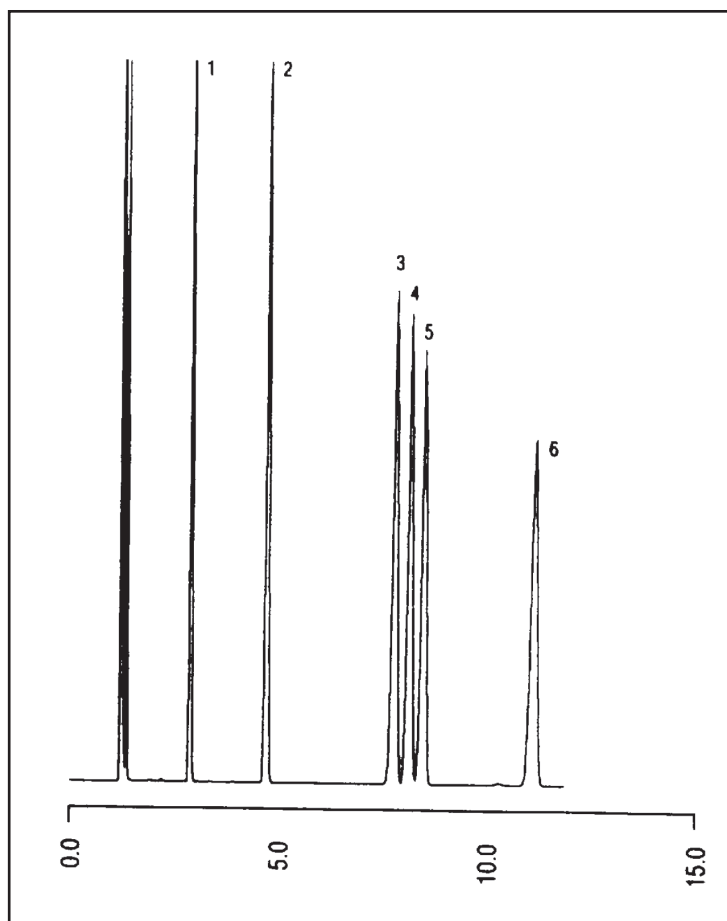

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: [tech@quadrexcorp.com](mailto:tech@quadrexcorp.com)