


# GC CAPILLARY COLUMN APPLICATIONS: ULTRA ALLOY™ STAINLESS STEEL

## Butter

COLUMN: UAC-TRG, Triglyceride Phase  
30M. x 0.25mm I.D. x 0.15µm film  
Cat. No.: UAC-TRG-30-0.15F  
Temperature: 250° (4°/min.) - 360°C

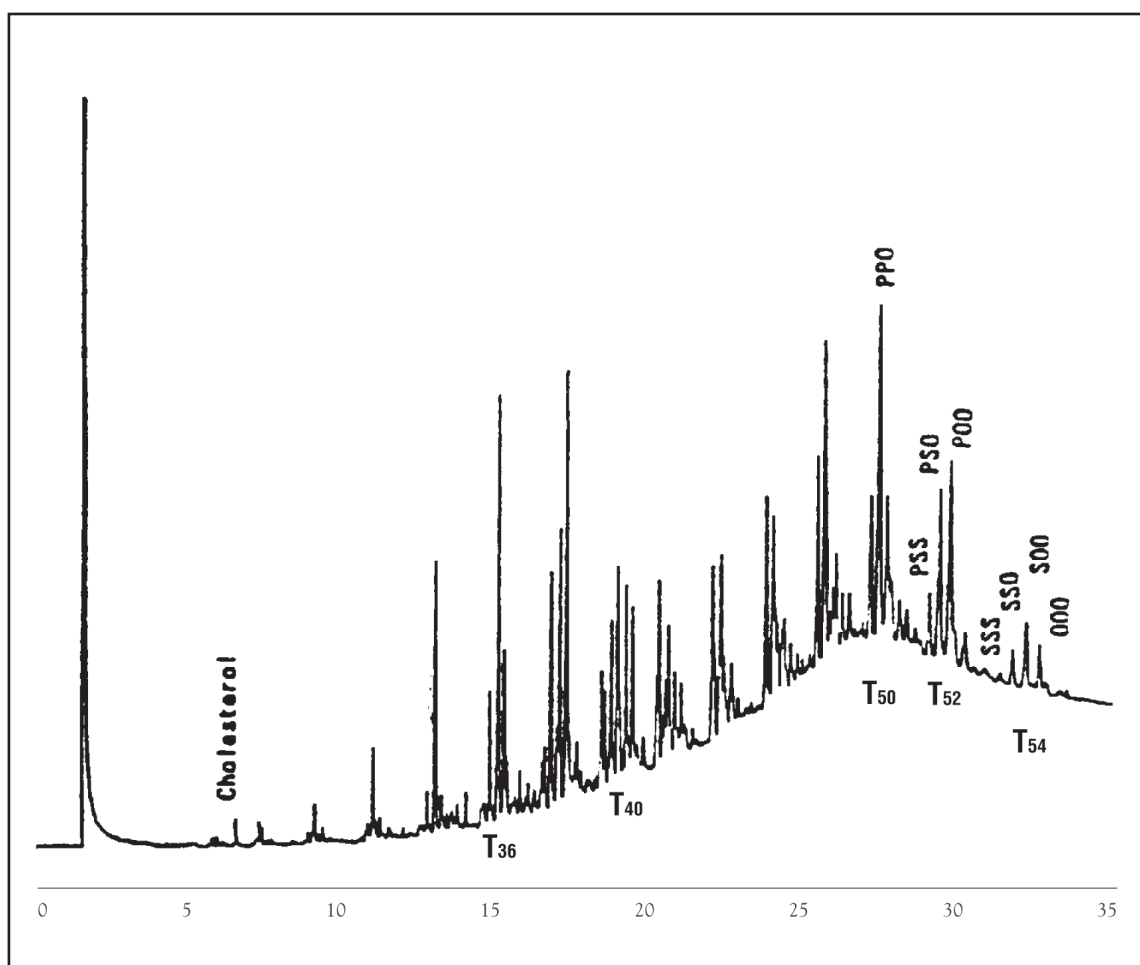

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: [tech@quadrexcorp.com](mailto:tech@quadrexcorp.com)