


# GC CAPILLARY COLUMN APPLICATIONS: ULTRA ALLOY™ STAINLESS STEEL

## C12 to C60 Hydrocarbons

COLUMN: UAC-DX30, Carborane siloxane  
15M. x 0.25mm I.D. x 0.15µm film  
Cat. No.: UAC-DX30-15-0.15F  
Temperature: 100° (20°/min.) - 450°  
Injector: 350°C  
Detector: 420°C, FID  
Carrier Gas: 30 cm/sec., Helium

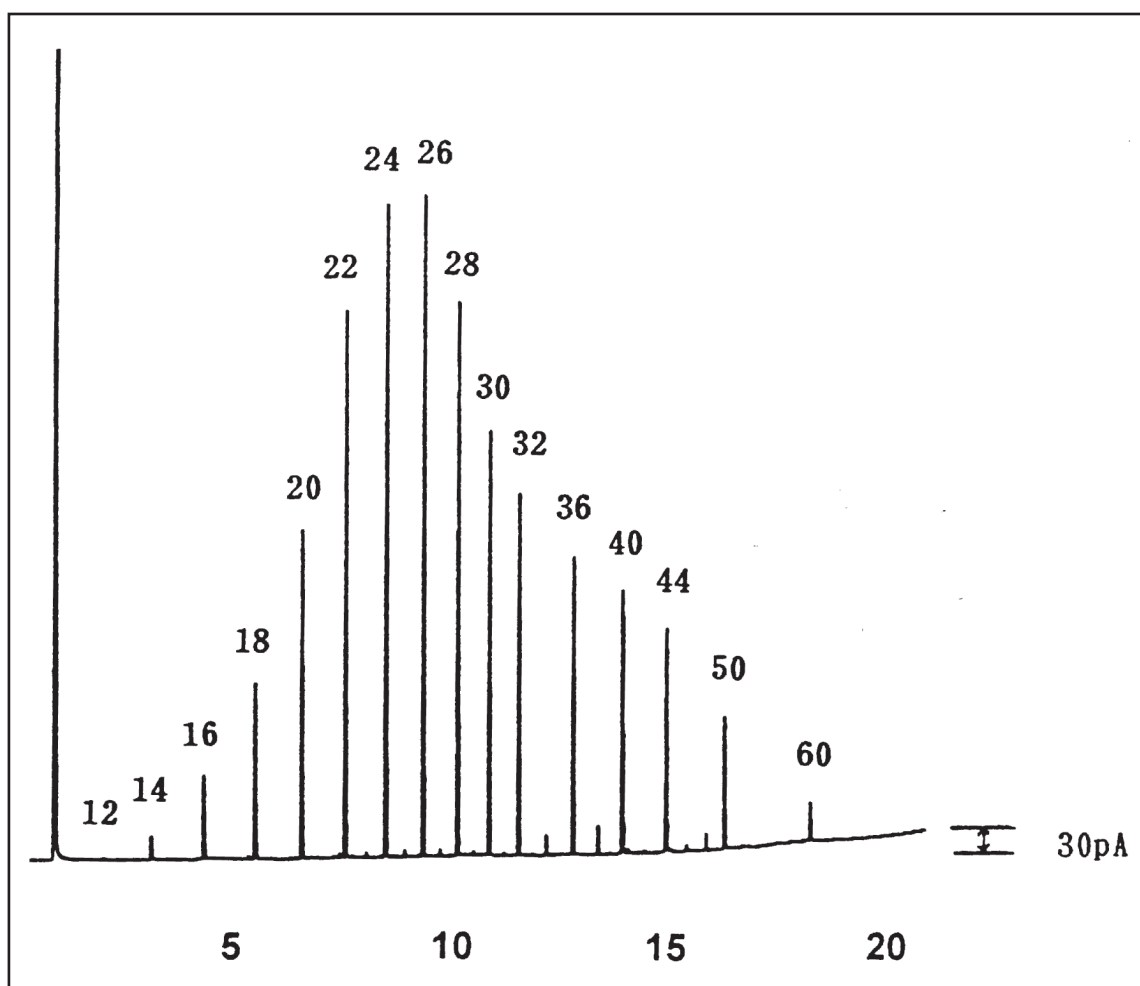

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: [tech@quadrexcorp.com](mailto:tech@quadrexcorp.com)