


# GC CAPILLARY COLUMN APPLICATIONS: ULTRA ALLOY™ STAINLESS STEEL

## Pesticides and Herbicides

COLUMN: UAC-1, Dimethylpolysiloxane  
15M. x 0.25mm I.D. x 0.25µm film  
Cat. No.: UAC-1-15-0.25F  
Temperature: 50° (1 min. hold) (20°/ min) -150° -  
(5°/min.) - 300°C  
Injector: 250°C  
Detector: 300°C, FID  
Carrier Gas: 0.8 ml/min., Helium

1. Chloneb
2. Pencycuron
3. Benfluralin
4. Simazine
5. Chlorothalonil
6. Propyzamide
7. Diazinon
8. Tolclofos-methyl
9. Terbutcarb
10. Fenitrothion
11. Chlorpyrifos
12. Captan
13. Pendimethalin
14. Methyl dymron
15. Isfenphos
16. Isoprothiolane
17. Napropamide
18. Butamifos
19. Flutolanil
20. Isoxathion
21. Mepronil
22. Pyridaphenthion
23. Iprodione
24. Bensulide

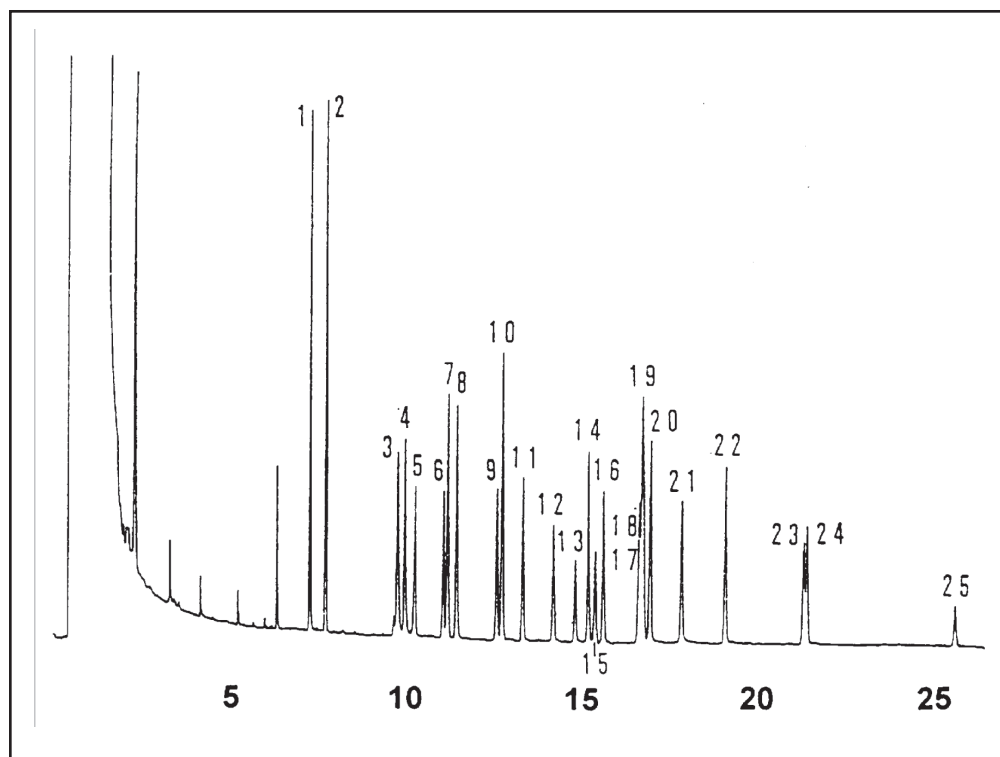

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: [tech@quadrexcorp.com](mailto:tech@quadrexcorp.com)