


# GC CAPILLARY COLUMN APPLICATIONS: ULTRA ALLOY™ STAINLESS STEEL

## Residue Solvents in Paint

COLUMN: UAC-CW, Carbowax  
30M. x 0.50mm I.D. x 1.0µm film  
Cat. No.: UAC-CW-30V-1.0F  
Temperature: 50° (3 min. hold)  
(5°/ min.) - 100°C  
Injector: 200°C  
Detector: 200°C, FID  
Carrier Gas: 3.4 ml/min., Helium

1. THF
2. ethylacetate
3. MEK
4. iso-C<sub>3</sub>H<sub>7</sub>OH
5. toluene

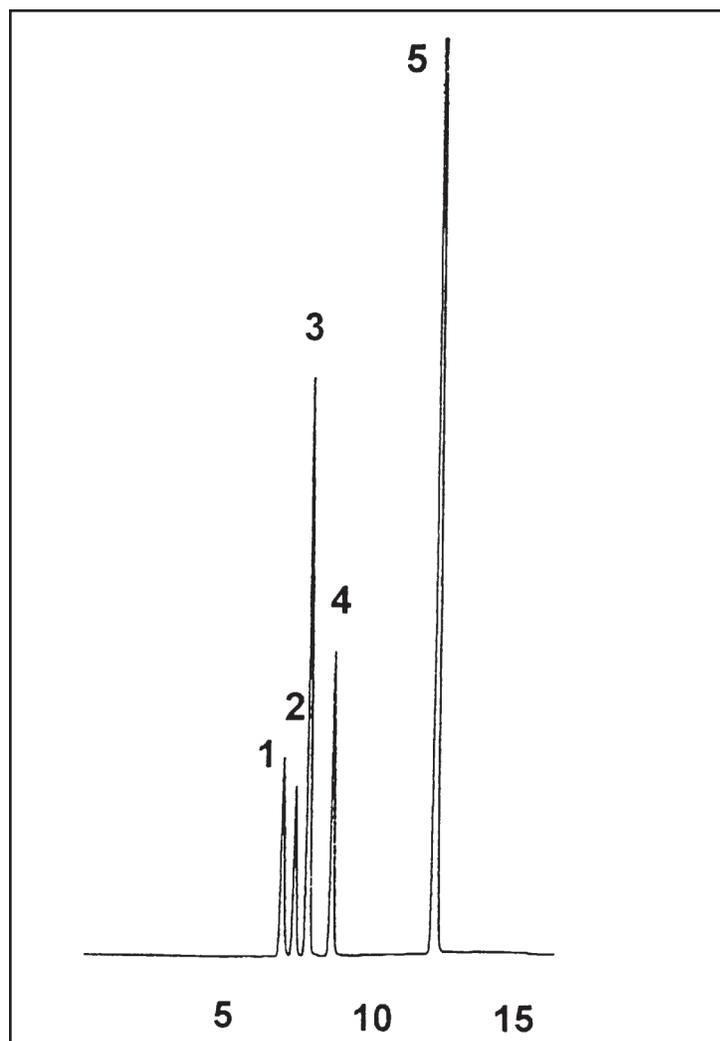

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: [tech@quadrexcorp.com](mailto:tech@quadrexcorp.com)